

# Metis™ Green

## Coronary Guide Wire

Any day, any case, any fingers



### Product Information

#### Features


- CoreToTip design
- Core/Coil technology
- 3 tip-load levels
- Hydroskin™ coating
- Straight / J-Tip
- Platinum/Nitinol spring coil
- Extendable

#### Sizing

- 0.014"
- Tip load 0.6, 1.0, 2.0 g
- 185 cm / 300 cm\*
- 3 cm platinum coil

#### Benefits

- Balanced flexibility
- Enhanced push
- Highly steerable
- Shapeable, re-shapeable
- Optimal tactile feedback
- Enhanced opacity
- Moderate support
- DOC wire compatible

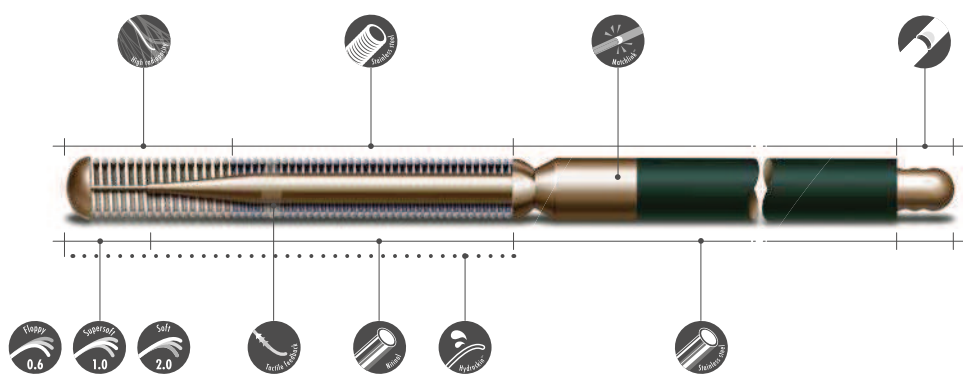

- Flattened nitinol core wire  
3 tip-load levels
- 3 cm opaque coil  
platinum/nitinol alloy
- 22 cm olucent coil  
stainless steel
- ACT™ Welding; 4 cm from tip,  
enhancing tactile feedback
- 25 cm Hydroskin™  
hydrophilic coating
- Nitinol  
core wire
- Matchlink™ Welding  
nitinol/stainless steel
- 135 cm stainless steel  
PTFE coated core wire
- 2,5 cm extension ribbon  
for DOC wire

# Ordering information

## Metis™ Green

METIS™ GREEN CORONARY GUIDE WIRE - 185 CM - 5 UNITS PER BOX

Product	Diam	Length	Design	Tip	Tip Load g	EXT*	Catalogue #	
Metis™ Green	0.014"	185 cm	Core To Tip	STR	Floppy	0.6	✓	CGW185C-0S
Metis™ Green	0.014"	185 cm	Core To Tip	STR	Supersoft	1.0	✓	CGW185C-1S
Metis™ Green	0.014"	185 cm	Core To Tip	STR	Soft	2.0	✓	CGW185C-2S
Metis™ Green	0.014"	185 cm	Core To Tip	J	Floppy	0.6	✓	CGW185C-0J
Metis™ Green	0.014"	185 cm	Core To Tip	J	Supersoft	1.0	✓	CGW185C-1J
Metis™ Green	0.014"	185 cm	Core To Tip	J	Soft	2.0	✓	CGW185C-2J

\*EXT : extension ribbon

- Additional Metis™ 300 cm can be ordered on demand at no additional cost. Minimum quantity of 6 units.

\* Extension Wire: Metis™ CGW (185 cm) can be extended by using compatible extension wire on the market. Please consult your local Comed Alliance Group member.

15 mm flattened Nitinol core wire  
shapeable and re-shapeable

