


# Metis™ Blue

## Coronary Guide Wire

Any day, any case, any fingers

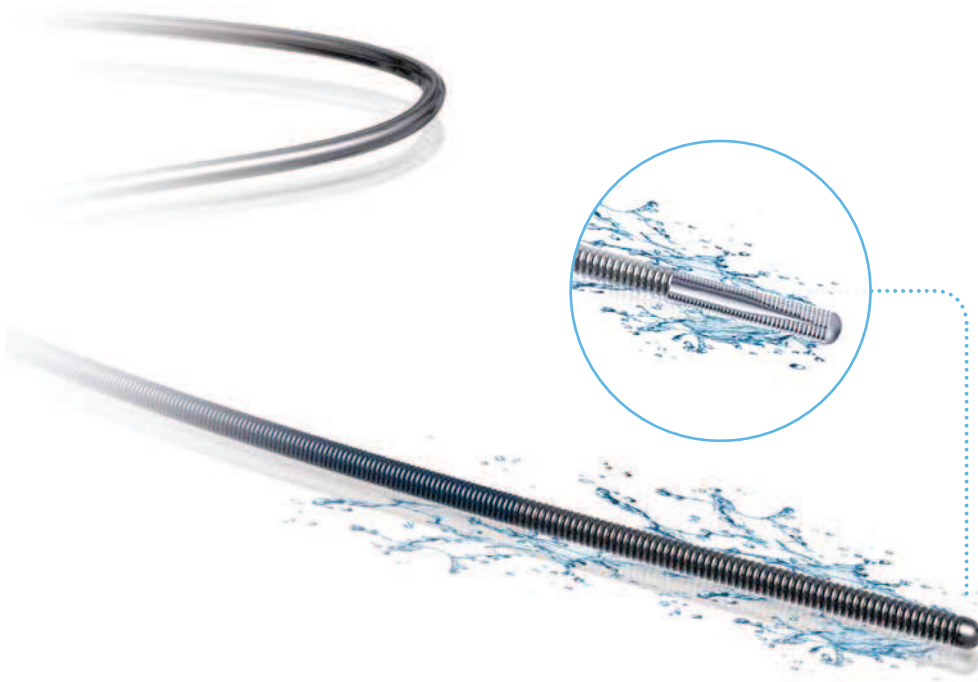

### Product Information

#### Features


- Shaping Ribbon design
- Core/Coil technology
- 3 tip load levels
- Hydroskin™ coating
- Straight / J-Tip
- Platinum/Nitinol spring coil
- Extendable

#### Sizing

- 0.014"
- Tip load 0.6, 1.0, 2.0 g
- 185 cm / 300 cm\*
- 3 cm platinum coil

#### Benefits

- High flexibility
- Optimal push
- Balanced steerable
- Shapeable, re-shapeable
- Optimal tactile feedback
- Enhanced opacity
- Moderate support
- DOC wire compatible

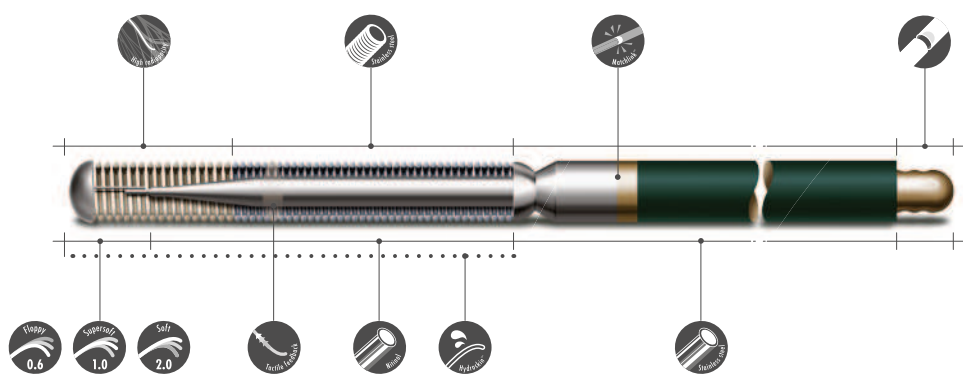

- Stainless steel shaping ribbon  
3 tip-load levels
- 3 cm opaque coil  
platinum/nitinol alloy
- 22 cm olucent coil  
stainless steel
- ACT™ Welding; 4 cm from tip,  
enhancing tactile feedback
- 25 cm Hydroskin™  
hydrophilic coating
- Nitinol  
core wire
- Matchlink™ Welding  
nitinol/stainless steel
- 135 cm stainless steel  
PTFE coated core wire
- 2,5 cm extension ribbon  
for DOC wire

# Ordering information

## Metis™ Blue Coronary Guide Wire

METIS™ BLUE CORONARY GUIDE WIRE - 185 CM - 5 UNITS PER BOX

Product	Diam	Length	Design	Tip	Tip Load g	EXT*	Catalogue #	
Metis™ Blue	0.014"	185 cm	Shaping Ribbon	STR	Floppy	0.6	✓	CGW185S-0S
Metis™ Blue	0.014"	185 cm	Shaping Ribbon	STR	Supersoft	1.0	✓	CGW185S-1S
Metis™ Blue	0.014"	185 cm	Shaping Ribbon	STR	Soft	2.0	✓	CGW185S-2S
Metis™ Blue	0.014"	185 cm	Shaping Ribbon	J	Floppy	0.6	✓	CGW185S-0J
Metis™ Blue	0.014"	185 cm	Shaping Ribbon	J	Supersoft	1.0	✓	CGW185S-1J
Metis™ Blue	0.014"	185 cm	Shaping Ribbon	J	Soft	2.0	✓	CGW185S-2J

\*EXT : extension ribbon

## Metis™ Blue Coronary Guide Wire

METIS™ BLUE CORONARY GUIDE WIRE - 300 CM - 3 UNITS PER BOX

Product	Diam	Length	Design	Tip	Tip Load g	Catalogue #	
Metis™ Blue	0.014"	300 cm	Shaping Ribbon	STR	Floppy	0.6	CGW300S-0S

- Additional Metis™ 300 cm can be ordered on demand at no additional cost. Minimum quantity of 6 units.

\* Extension Wire: Metis™ CGW (185 cm) can be extended by using compatible extension wire on the market. Please consult your local Comed Alliance Group member.

13 mm stainless steel  
flattened Shaping Ribbon  
shapeable and re-shapeable

